Specification

Capacity:3kg-180kgs
 Division:1D=100g/0.2lb
 Unit:KG/LB/ST
 Battery:2*AAA battery(not include)
 Best working temperature:5~35℃
 Storage temperature:-10~60℃
 Moisture Environment:≤90%RH

Notices

KG/LB Switch button: There is a button on the back of the weight scale.
 1. Please stand on the scale in the right direction.
 2. Please stand on the scale no less than 10 seconds if you measure weight.
 3. Please confirm gender, height, weight data is correct

Warnings

1. Please keep the scale away from wet floor, don't use the scale if the glass is wet in case of slip.

2. Please stand on the middle of the scale in case of turnover.
 3. The platform of the scale is made of glass, please avoid sharp or harsh things scratching the surface, excessive forces or sudden shocks, unexpected fallen-off will break the glass and injure the user.
 4. The scale is not waterproof. Keep the surrounding environment dry, do not submerge it in water.
 5. Use soft cloth with neutral cleaning liquids to clean the surface, keep away contacting with organic solvent nor chemicals.
 6. When the batteries are in low power status, please remove the batteries in case of leaking liquids corrode the scale.

English

Step 1 Install App

Use your device to scan QR Code or search OKOK from App Store/Google Play to install App.



- 1/ Switch on wireless
- 2/ Phone or Tablet connects to Net (WiFi, 3G or 4G)

Step 3 Weighing

Enter into "Dynamic Weight" interface

Stand on the scale and synchronize weight data onto APP

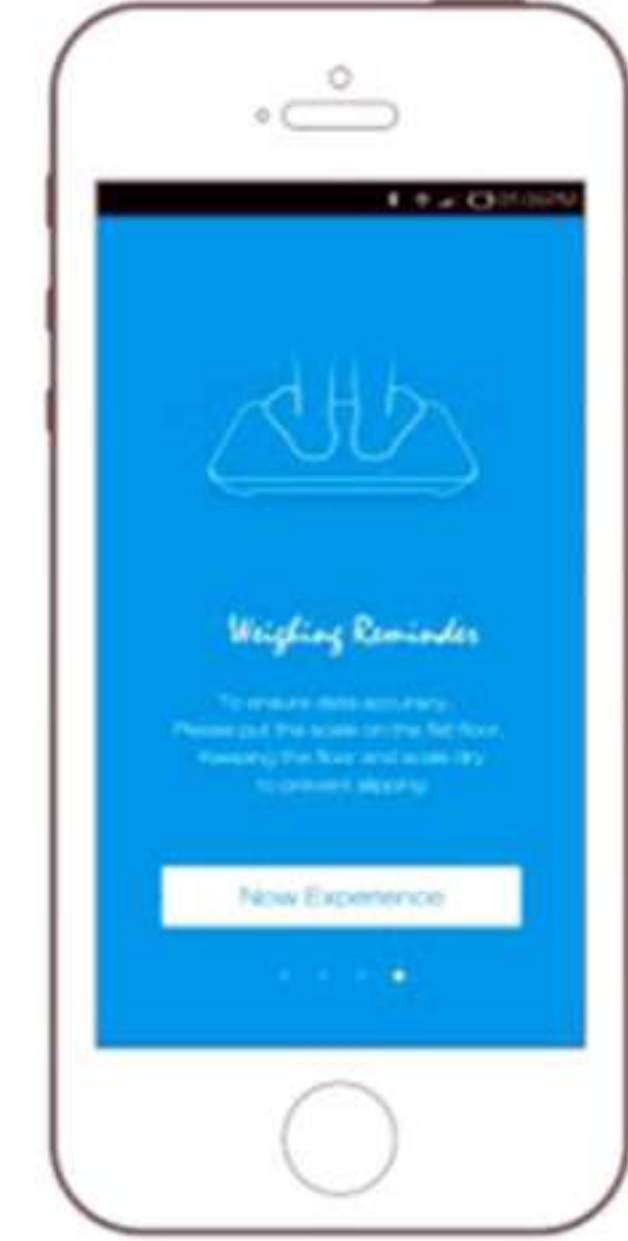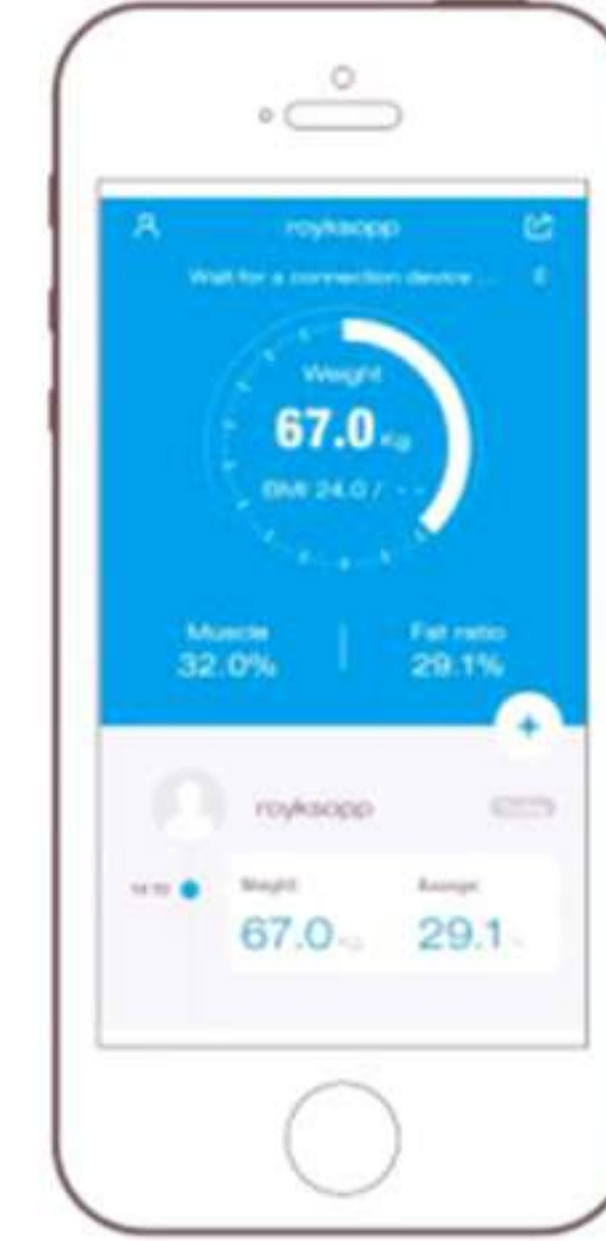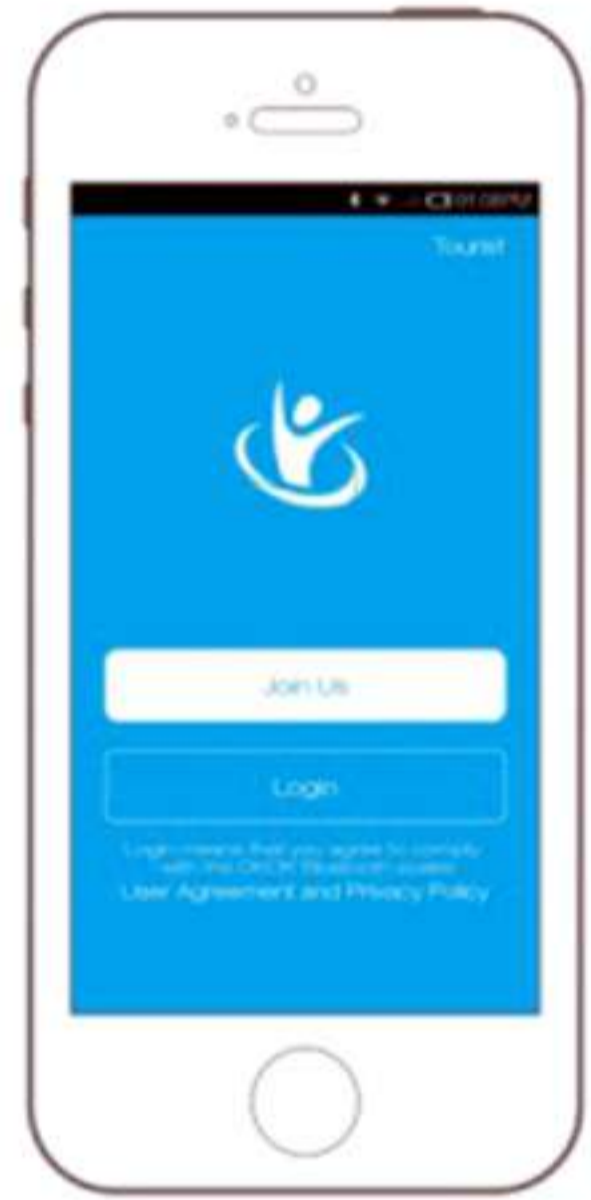
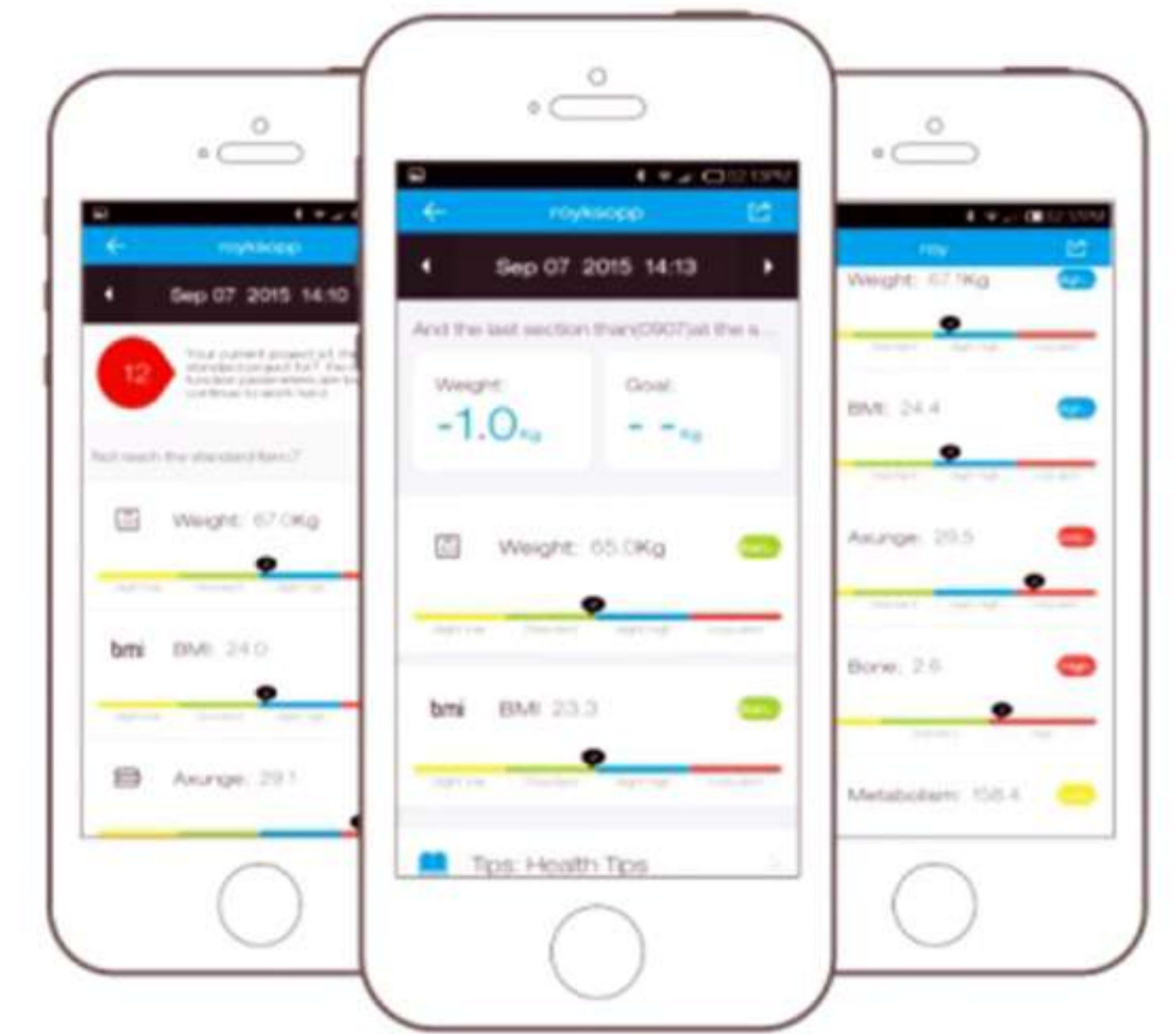
Step 2 Pair the devices, register and complete user information

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